Rotadent ProCare®

The Professional Rotary Toothbrush

Instructions for Use

English • Español • Français • Italiano •
Português do Brasil • Deutsch • 日本語 • العربية

IMPORTANT

- Charge your Rotadent for at least 24 hours prior to initial use.
- Use a gel-based toothpaste to best preserve the MicroFilaments™.

**To order Rotadent
brush heads or accessories:**
Please contact your Dental Professional

TABLE OF CONTENTS – TRANSLATIONS

English	2
Español	10
Français	20
Italiano	30
Português do Brasil	40
Deutsch	50
日本語	60
العربية	70

TABLE OF CONTENTS

IMPORTANT SAFEGUARDS	2
DANGER	2
WARNINGS	2
MEDICAL WARNINGS	4
BATTERY RECYCLING	5
 Rotadent ProCare® Attachments	 5
Brush Heads	5
MicroFilaments™	5
Color Rings	5
 Preparing for Use	 6
Charging Your Rotadent ProCare.....	6
Assembling Your Rotadent ProCare	6
 Brushing with Your Rotadent ProCare	 7
Brushing Modes	7
Changing Speeds	7
QuadTimer	7
LED Indicators	7
 Brushing Technique	 8
Tingling Sensation	8
Tips for Special Oral Care Needs	8
 Cleaning Your Rotadent ProCare	 9
 Storing Your Rotadent ProCare	 9
 Troubleshooting and Guidelines	 9

Rotadent ProCare®

The Professional Rotary Toothbrush

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

PLEASE READ ALL INSTRUCTIONS BEFORE USING

DANGER

To reduce the risk of electrocution:

1. Do not use this appliance while bathing.
2. Do not place or store the charging base where it can fall or be pulled into water or other liquid.
3. Do not reach for the charging base if it has fallen into water or other liquid. If the charging base or cord falls into water or other liquid, unplug it immediately.
4. Never use a charging base that has a damaged cord or wall plug. If damaged, please contact your distributor for a replacement charging base.

WARNINGS

To reduce the risk of burns, electrocution, fire, or injury to persons:

1. Use this appliance only for its intended use as described in this manual or as recommended by your Dental Professional.

2. This appliance is not to be used or cleaned by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
3. Children shall not play with this appliance.
4. Do not use attachments other than those provided or recommended by DenMat.
5. The power cord connected to the charging base cannot be replaced. If the power cord is damaged, please dispose of the charging base and contact distributor for a replacement.
6. Do not force the wall plug into an outlet.
7. Keep the cord away from heated surfaces.
8. Never operate this appliance while drowsy.
9. This appliance is for household use only.
10. This appliance contains no user-serviceable parts. Refer to the **Troubleshooting Guidelines** section if the appliance no longer works properly.
11. Recharge the battery only with the charging base that has been supplied with the appliance.

12. If you are traveling, ensure that the voltage indicated on the charging base corresponds to the local voltage before connecting the appliance.

Charging Base Item #	Electrical Rating
94000CB	120 V, 50-60 Hz, 150 mA, CA uniquement
94000CBUK 94000CBEU	230V, 50-60 Hz, 150 mA, CA uniquement

MEDICAL WARNINGS

1. Consult your Dental Professional before you use this appliance.
2. Contact your Dental Professional if prolonged discomfort is experienced during use.
3. This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician prior to use.

SAVE THESE INSTRUCTIONS

4. Consult your physician if you have other medical concerns surrounding the use of this product.

Battery Recycling

When the handle has reached the end of its useful life, please obey local recycling laws by recycling the rechargeable battery. Do not dispose of the battery with normal household waste.

Rotadent ProCare® Attachments

Brush Heads

DenMat offers three types of brush heads for the Rotadent ProCare. It is recommended that you replace your brush head every three months to achieve optimal results.

- **Hollow-Tip**

The hollow-tip brush head is designed for cleaning and polishing surfaces of the teeth and sweeping under the gumline. Its flat, cup-like design makes it ideal for absorbing and applying oral medication, as prescribed by your Dental Professional.

- **Short-Tip**

The short-tip brush head is the most universal brush head. Like the hollow-tip, it is effective for cleaning and polishing surfaces of the teeth and is also able to deliver oral medication more deeply into periodontal pockets and between teeth.

- **Long-Tip**

The long-tip brush head is designed for patients with large spaces between teeth or for patients with orthodontics or prosthetics. It is most effective for targeting hard-to-reach areas and sweeping plaque and debris away from braces, multi-unit fixed bridges or between teeth.

MicroFilaments™

Each Rotadent ProCare brush head features 4,577 MicroFilaments™ that reach below the gumline and into tight interproximal spaces. Because these MicroFilaments are so thin and dense, you may find that traditional toothpaste is too heavy and thick to use with your Rotadent ProCare.

IMPORTANT: Use of a gel-based toothpaste is recommended in order to maximize the useful life of the MicroFilaments™.

Color Rings

Each of the three brush heads comes with a color ring attached to it. These color rings are interchangeable and can easily be removed and replaced. In households where one Rotadent ProCare is being shared amongst several users, these color rings can be used as a way of keeping track of which brush head belongs to each user.

Preparing for Use

Charging Your Rotadent ProCare®

The Rotadent ProCare is powered by a rechargeable nickel-metal hydride (NiMH) battery. When fully charged, the Rotadent ProCare may last up to two weeks (when used for two minutes, twice a day), depending on brushing conditions.

Initial Charging

- Plug the power cord of the charging base into a working outlet.
- Place the Rotadent ProCare into the charging base. The charging indicators will blink.

IMPORTANT: In order to achieve a full charge, charge your Rotadent ProCare for at least 24 hours prior to initial use.

Recharging

When the power level of your Rotadent ProCare noticeably decreases, it is time to recharge the battery. Follow the steps listed under the **Initial Charging** section.

Once the battery has reached a full charge, the LED indicators will stop blinking (as shown in the **LED Indicators** section on the next page). At that point, there is no risk of overcharging the battery by leaving the handle on the charging base. For everyday use, you can safely leave the handle on the charging base with the charging base plugged into a wall outlet, without affecting the useful life of the battery.

IMPORTANT: The handle is sealed for your safety and the battery is not accessible.

Assembling Your Rotadent ProCare

Color Ring Attachment

- Orient the color ring so that the taper matches the taper of the neck.
- The color ring will be oriented on the neck so that the internal detail matches the external detail on the lower portion of the neck.
- Attach the color ring by hooking one edge on the bottom of the brush head and stretching it into place.
- The color ring can alternately be attached by sliding it over the top of the brush head and down the neck, into position.

Brush Head Attachment


- Align the brush head so that the MicroFilaments™ are facing outward from the front of the handle (the same direction as the on/off button).
- The inner components of the brush head should align easily with the brush head connection at the top of the handle.
- Push the brush head down onto the handle until it locks into position. You should hear a very soft “click” sound, confirming that the brush head is securely in place.

Brushing with Your Rotadent ProCare®

Brushing Modes

The Rotadent ProCare offers two different brushing modes to suit your comfort level:

- **Normal** — Recommended mode for superior teeth cleaning.
- **Slow** — Provides gentle and thorough cleaning for sensitive gums and teeth at a reduced speed.

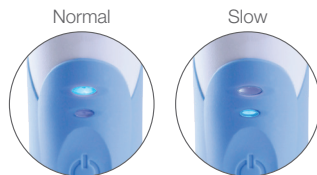Changing Speeds

- When pressing the on/off button, your brush automatically starts in Normal mode, which is indicated by the larger LED light.
- To switch to Slow mode, press the on/off button again until the smaller LED light becomes illuminated.
- To turn your brush off, press the on/off button. If the Rotadent ProCare has been running at either speed for more than 30 seconds, the next press of the on/off button will turn the unit off. If it has been running for less than 30 seconds, each press of the on/off button will change to the next mode selection. Pressing the on/off button while in the Slow mode will turn the unit off.
- The Rotadent ProCare may also be turned off by placing the brush in the charging base.

QuadTimer

- The Rotadent ProCare features the two-minute QuadTimer to assist you in reaching the recommended two-minutes of brushing.
- The QuadTimer helps you brush evenly and thoroughly with a pause at 30, 60, 90 and 120 second intervals.
- The brush will run for an additional 30 seconds after the fourth pause to allow you to return to areas that may have been missed.
- After 2 minutes and 30 seconds (150 seconds), the brush will shut off, indicating that you have reached the full brushing cycle.

LED Indicators

Brushing Modes		Immediately After Use		Handle is on Charging Base and Plugged In	
Normal	Slow	Fully Charged	50% or Less	Battery is Charging	Battery is Fully Charged


IMPORTANT: After pressing the on/off button to turn the unit off, the LED charge indicator will remain illuminated for a period of 30 seconds, indicating the level of charge remaining.

Brushing Technique

Rotadent ProCare® requires a tooth-by-tooth brushing technique. The 360-degree rotation may cause the brush head to slip off of the teeth during use, particularly when first getting used to this technique. To help correct this, stabilize the handle with a firm grip.

Next, follow the margins of the gumline to ensure proper cleaning of all surfaces, as described below:

- Hold your Rotadent ProCare so the brush head meets your gumline and the tooth surface at a 90-degree angle.
- To thoroughly clean upper teeth, make “M” motions to massage the gumline.
- To thoroughly clean lower teeth, make “W” motions to massage the gumline.
- To make sure you brush evenly, divide your mouth into four sections (or quadrants). Brush each quadrant thoroughly for 30 seconds, as indicated by the QuadTimer.

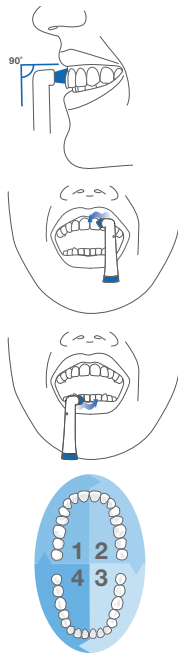TECHNIQUE TIP: To most-effectively clean in between teeth or around braces, utilize the “M” and “W” motions described above, using the long-tip brush head.

Tingling Sensation

Some individuals may experience a mild tingling sensation when first using the Rotadent ProCare. This should stop after a few days of use.

Tips for Special Oral Care Needs

Braces

With braces, it is not always possible to hold your Rotadent ProCare at a 90 degree angle. If this is the case, be sure to clean the wires and brackets adequately as plaque builds on those areas as well as on teeth and gums.

Dental Restorations

The Rotadent ProCare is safe for use on dental restorations (such as veneers, crowns, implants, fillings and bridges). Use the Rotadent ProCare® as you would on all of your other normal tooth surfaces, making sure to work the brush tip around all surfaces of your dental restorations.

Cleaning Your Rotadent ProCare®

Cleaning Your Brush Head

To clean the brush head, remove the brush head from the handle and rinse separately under warm running water.

Cleaning Your Handle

Use a damp cloth to wipe the entire surface of the handle. Make sure the handle is completely dry before placing the handle back on the charging base.

Cleaning Your Charging Base

Before cleaning your charging base, make sure to unplug the power cord. Place the charging base on a safe and dry surface. Only use a dry cloth or damp cloth to wipe the surface of the charging base. Do not plug the charging base into the wall outlet until it is completely dry.

IMPORTANT: Do not clean the charging base by immersing it under water.

Storing Your Rotadent ProCare

If your Rotadent ProCare® will not be used for an extended period of time, unplug the charging base, clean the appliance as described in the **Cleaning Your Rotadent ProCare** section, and store it in a cool, dry place away from direct sunlight.

Troubleshooting Guidelines

PROBLEM	SOLUTION
The handle is not working...	a) The battery may need to be charged. Recharge the handle for at least 24 hours. b) The outlet may not be live. Try a different power outlet. The electricity supply to outlets in bathrooms may be cut off when the light is switched off. c) The handle is on the charging base. Remove it from the charging base (for safety reasons, the toothbrush will not work in the charging base). d) Brush head may be damaged. Replace with a new brush head.
The handle is running slowly...	a) Battery charge may be low. Put the handle into the charging base to recharge. b) It may be time to change your brush head. Replace with a new brush head.
The brush head is noisy...	a) The brush head may not be locked into position. Check the brush head installation per the instructions in this manual. b) Excessive toothpaste buildup could be keeping the brush head from working correctly. Clean the handle and brush head, and then try to lock the brush head in place again. c) Replace with a new brush head.

If you are still experiencing technical issues with your Rotadent ProCare, please contact your local DenMat distributor for further assistance.

Rotadent ProCare®

The Professional Rotary Toothbrush

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©2018 Den-Mat Holdings, LLC. All rights reserved. Reservados todos los derechos. Tous droits réservés. Tutti i diritti riservati. Todos os direitos reservados. Alle Rechte vorbehalten.
Den-Mat Holdings, LLC